


AIR COMPRESSION THERAPY SOLUTIONS


All human body is supplied by peripheral circulation of blood and lymph.

Our systems help patients to recover from their disease by improving blood and lymph circulation and metabolism.

Our systems consist of air chambers and an air pump offering air compression massage to legs, arms or waist.

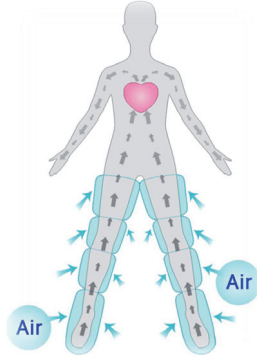


• **28430 MAYLEA GIMA DIGITAL COMPRESSION MASSAGE with 2 legs (98x38 cm) - 4 chambers, 1 single and 1 double hose and 2 plastic feet guides**
Skip function to choose directly the parts to be massaged. Fully digital management. Autocalibration functions.
Manual: GB, FR, IT, ES. This system can operate with **3 compression modes (A, B, C).**

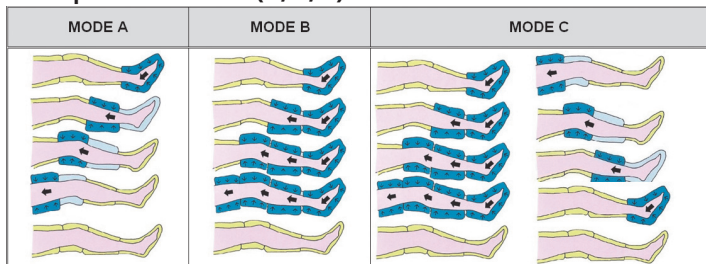

TECHNICAL SPECIFICATIONS	
- pressure sensor	- one touch connection hose plug
- digital timer (10, 20, 30, 40 min)	- body size: 190x250x210 mm
- digital pressure gauge (10-180 mmHg)	- body weight: 3 kg

- **28435 LEG 93x32 cm**
- **28434 LEG EXTENSION ZIPPER - 10 cm extra extension**
- **28440 DOUBLE HOSE - to use two leg cuffs at the same time - spare**
- **28436 ARM CUFF unique size - 86x29 cm**
- **28437 WAIST CUFF 130x38 cm**

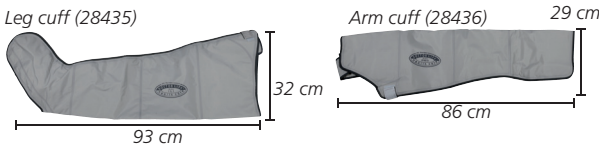


• **28441 MK400L PROFESSIONAL COMPRESSION SYSTEM with 2 legs (XL), and 2 single hoses**
MK400L is a professional air compression therapy system with essential functions for different kinds of patients. This system provides 3 compression modes (wave, squeezing, double squeezing), reverse compression settings, skip functions to select directly part to be pressurized and compression speed controller. Touch key pad display. Manual: GB, IT, ES.

TECHNICAL SPECIFICATIONS	
- Pressure range: 10-200 mmHg	- Use time: 1-90 minutes
- Automatic pressure detection	- Body size: 415x310x160 mm
- Operating voltage: 220-240 V, 50/60 Hz, 50 VA	- Body weight: 6.5 kg



- Medical indications**
- DVT (Deep Vein Thrombosis)
 - Lymph edema
 - Post-paralytic
 - Limb convulsion
 - Venosclerosis
 - Chronic venous insufficiency (CVI)
 - Post-mastectomy
 - Post-trauma
 - Spinal cord injury
 - Varicose veins
 - Sciatica
 - Sprains
 - Diabetes
 - Arterial insufficiency

- Home care indications**
- Lymph edema
 - Diseases caused by interrupted blood circulation
 - Cellulitis fat dissolution
 - Pregnant woman's limb treatment
 - Numb feet or hands
 - Cold feet or hands
 - Settlement of stress
 - Exercise effects
 - Insomnia
 - Healthier body and elastic skin
 - Old age, weakness and fatigue

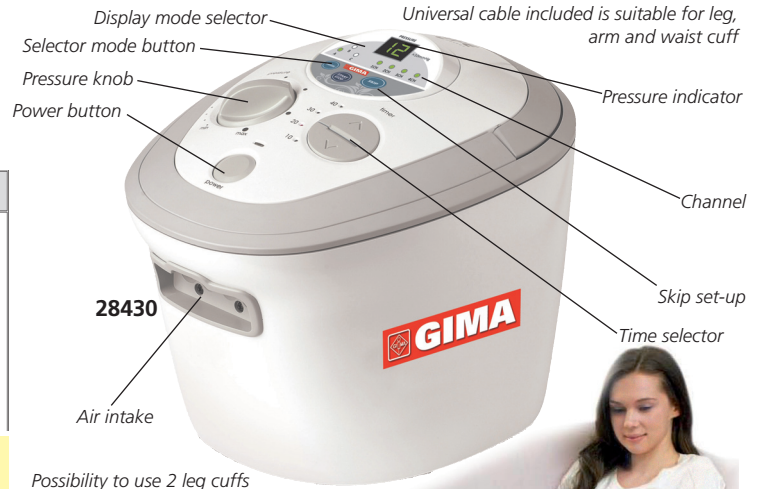

Possibility to use 2 leg cuffs at the same time with double hose (code 28440)



PROFESSIONAL AIR COMPRESSION SYSTEM




- **28442 LEG CUFF L - 6 CHAMBERS - spare**
- **28443 LEG CUFF XL - 6 CHAMBERS - spare**
- **28444 LEG EXTENSION ZIPPER XL - 10 cm extra extension for 28443**
- **28445 SINGLE HOSE - spare**
- **28446 ARM CUFF - 6 CHAMBERS - universal**
- **28447 WAIST CUFF - 6 CHAMBERS - min-max girth 94-125 cm**

A mode (Wave)
The selected air chamber inflates and deflates by one chamber sequentially from foot up to thigh.

B mode (Squeezing)
After all the selected air chambers inflate sequentially from foot up to thigh the inflated chambers deflate at once.

C mode (B mode two times + Reverse A mode 1 time)
The selected air chambers inflate and deflate by two chambers sequentially from foot up to thigh.
Reverse setting: it allows to change automatically the selected mode to reverse mode.

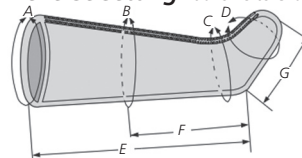

Cuff Size	Leg dimension (cm)						
	A	B	C	D	E	F	G
L	64,0	49,5	40,5	34,5	75,0	38,0	28,0
XL	73,5	55,0	42,0	35,5	85,0	43,0	31,0