

# JD-2600A AUTO LENSMETER

NEW



- Adjustable LCD display screen

Tilt color screen provides you with clearview from any position such as sitting and standing.

- Sensitive touch LCD display screen

Just touch the screen can finish your measurement, your operation could be easier and more convenient.

This new generation auto lensmeter is based on Hartmann Sensor technology. Advanced simultaneous measurement of 108 points provides a faster measurement with greater accuracy and reliability.





### ● Measurement screen

Easy to use with recognisable icons. The JD-2600A provides a user-friendly graphical interface for an immediate understanding.

### ● Progressive measurement

Automatic measurement of progressive lenses. A graph helps the operator to improve accuracy and speed in the measurement of near and far vision.

### Range of measurement

|              |                               |
|--------------|-------------------------------|
| Sphere       | 0~±25D 0.01/0.12/0.25D step   |
| Cylinder     | 0~±9.99D 0.01/0.12/0.25D step |
| Axis         | 0°~180° (1°step)              |
| Add.         | 0~+9.99D 0.01/0.12/0.25D step |
| Prism degree | 0~15Δ 0.01/0.12/0.25Δ step    |

### Measurement mode

|                |                                           |
|----------------|-------------------------------------------|
| Cylinder       | +,+/-,-                                   |
| Prism          | X-Y,P-B                                   |
| Contact lens   | Soft/hard                                 |
| Measuring mode | Single/progressive/ automatic recognition |

### Specifications

|                      |                           |
|----------------------|---------------------------|
| Diameter of lens     | Φ20~108mm                 |
| PD                   | 40~90mm, 0.5mm step       |
| Speed of measurement | 0.1s                      |
| Display              | 7" TFT LCD 758×480 pixels |
| Printer              | Thermal Printer           |
| Dimension            | 222(D)×161(W)×372(H)mm    |
| Package              | 515(H)×320(W)×390(D)mm    |
| Weight               | ~5KG                      |
| Power supply         | 100~240V 50~60HZ          |

Specifications are tentative and subject to change

## JD-2600A AUTO LENSMETER

TAIYUAN XINYUAN HIGH-TECH CENTER NORTH UNIVERSITY OF CHINA

Tel: +86-351-3559336 | +86-351-3559527 Fax: +86-351-3921099

Http://www.xinyuan6000.com E-mail: fa6000@sina.com

Distributed By:


## Lab Medica Systems Pvt. Ltd.

705, Sakar 5, Nr. Mithakhali Railway Crossing,  
off Ashram Road, Ahmedabad-380 009, Gujarat, INDIA

Tel: +91-79-26574383 Fax: +91-79-26576828


M.: +91 9925017912, +91 9825024708, +91 9099926039

Email: info@labmedicasystems.com www.labmedicasystems.com

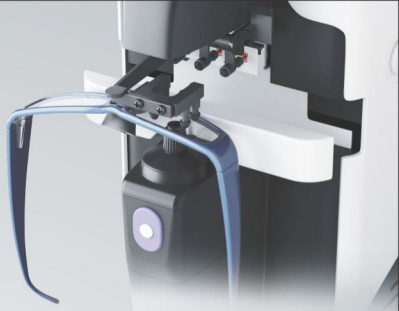

- Automatic recognition of single and progressive lens
- Hi-speed UV measurement

Just hold the lens through the sensor, the display will easily indicate the UV measurement by the progress bar.

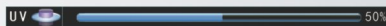

- Built-in thermal printer

Quick data print with built-in printer. Easy to reload printing paper.

- Built-in UV tester

