


# PROFESSIONAL SKINFOLD CALIPERS: MECHANICAL AND DIGITAL

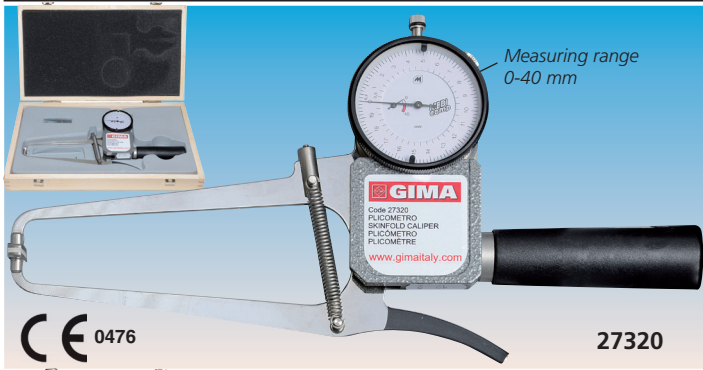

- 27320 PLICOMETRO - Professional skinfold caliper - mechanical
- 27346 PLICOMETRO - Professional skinfold caliper - digital with PC cable for data transfer USB

This high precision instrument measures the thickness of a fold of skin with its underlying layer of fat giving physicians the situation of the subcutaneous tissues of the patient.

### Possibility of re-calibration

Produced to meet International Standards, the Gima "Skinfold caliper" has springs which exert a pressure of 10 g/mm<sup>2</sup> on the skinfold and an accurate scale which measures thickness in mm. With possibility of re-calibration of the pressure it solves the major problem of traditional skinfold caliper that is reduction of the tension of the spring in the time.

### Easy to use and read

With double reading scale (code 27320) or digital display (code 27346), it is supplied with picture showing advised positions for measurements. Tables will help obtain correct % of body fat from measured data and make a correct evaluation of the patient.

### Stylish and user-friendly

A special nylon finish makes the device feel skin, avoiding trouble of contact with metal. Made in metal and sold in an attractive wooden box. Made in Italy, provided of medical device CE certification. User manual in English, French, Italian, Spanish.

### Data management on PC (only 27346)

For immediate data transfer on most programmes (Excel, Word, Open Office, ...). Requirements: Windows 2000, XP (Service Pack 2), 7. USB 2.0 minimum.



## BMI TAPE



• 27341 BMI - Tape Measure  
Pocket retractable tape measure 150 cm. BMI (Body Mass Index) calculator on front panel.

• 27342 BMI - Tape Measure  
Economical version of 27341. Available in 150 days for minimum orders of 1,000 pcs.



## FAT-1 SKINFOLD CALIPER

- 27344 FAT-1 SKINFOLD CALIPER


Economic but precise and accurate skinfold caliper for professional use, made in Italy. Measures the thickness of skin folds over the body, based on the characteristics developed by Tanner and Whitehouse:

- pressure on the fold: approx 10 g/mm<sup>2</sup>
- contact surface on the fold: approx 6x11 mm

It allows a correct evaluation of the nutritional status and the sectoral distribution of adipose tissue. Raise a fold of skin with your fingers and measure its thickness by the forceps: on the top edge a scale in millimeters and an arrow measure the skin fold (triceps, subscapular muscles, etc...) with an accuracy of about 1/2 mm.

The coloured bands placed near the scale allow a simultaneous comparison of the measured value with referential percentages related to men and women. This information is included in the manual provided with the unit and helps an easy interpretation of the values. Advised in dietetics, cardiology, sports medicine, aesthetic medicine, gynaecology, pediatrics, ecc...

Multilanguage manual: GB, FR, IT, DE, ES, GR.

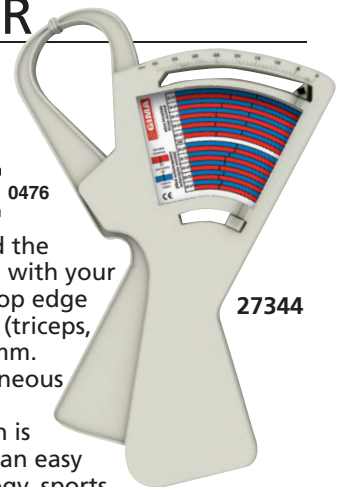

## BODY TAPE

- 27343 BODY TAPE MEASURE  
Track your body measurements and progress.

Measures any body part: arm, thigh, calf, chest, waist, hip....

Measures in inches and centimeters; total measurement length:

60 inches / 152 cm.

Stylish design, easy measurement.

Sturdy vinyl tape. Push-button retraction and locking feature ensure easy measurement and superior accuracy.



## PROFESSIONAL BODY - HEAD TAPE

- 27347 PROFESSIONAL BODY TAPE
- 27350 PROFESSIONAL BODY TAPE with pediatric gauge box

Patented, professional body tape measure approved by the Italian Ministry of Health for adult and pediatric use. Tape is coloured on both sides to help to identify metabolic risk in abdominal measure for male and female. Length of tape: 8-150 cm. Thanks to its rounded shape, it can be used also as a head measure tape. Packed in single box.

### Body tape with pediatric gauge box

Specific tape for pediatric use. Supplied in a box with a ruler to compare the measured head circumference with the standard percentiles of children's growth from 0 to 36 months. Suitable for both sexes.



Multilanguage manual: GB, FR, IT, DE, ES, GR, PT, Arabic.

